



Therapeutic Periodontal Supplement

One product for all periodontal issues

The Periodontitis




It is the **sixth most prevalent non communicable** chronic disease worldwide



60% of the population is susceptible to periodontal disease, with severe forms affecting **over 10% of the global population.**



It predominantly affects individuals in their **forties** with a low inclination towards dental hygiene.

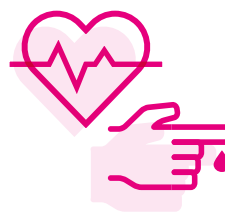

It is a systemic inflammation capable of **complicating cardiovascular diseases and diabetes.**



Our mission

Our priority is to contribute to prevention and ensure individuals affected by periodontal disease **have access to highly effective treatments with an optimal cost-benefit ratio.**

To achieve this, we have developed a **product to complement supra- and sub-gingival cleaning (SRP) as part of an advanced therapeutic program.**




WHAT IS CHLO-SITE ?

Chlo-Site is the innovative XANTHAN gel with BI-CHLORHEXIDINE designed to optimize periodontal care.

Chlo-Site has an occlusive and mucoadhesive action, remaining in the pocket by adhering to the tissues in a progressive 20+ day inhibition process

This feature plays a key role in preserving the hygiene established by the operator, amplifying the benefits of the treatment through the significant reduction of bacterial survival. This results in an optimal environment that effectively promotes the healing process.

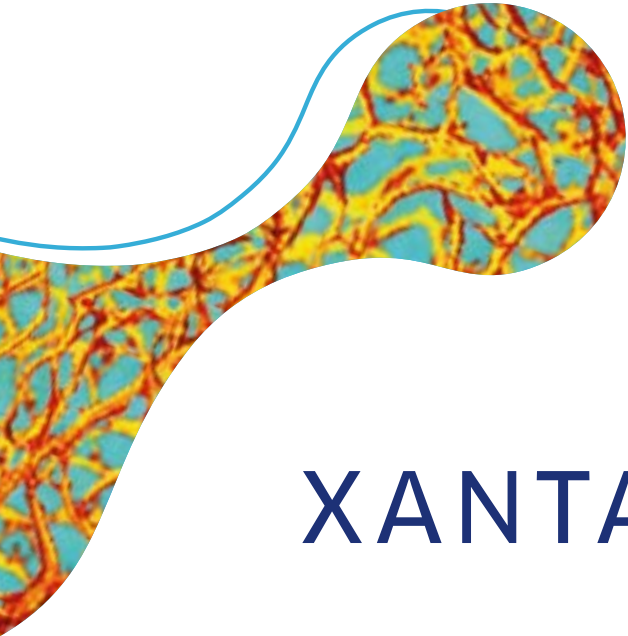

XANTAN + BI-CHLOREXIDINE

XANTAN is a polymer that forms a three-dimensional pseudo-plastic network, capable of suspending and retaining various substances, which are gradually released based on their physical and chemical characteristics. This makes it an excellent substrate for the formation of a stable gel and usable as a carrier for chlorhexidine.


CHLOREXIDINE is capable of preserving the application site from risks of bacterial contamination. In Chlo-Site, bi-chlorhexidine is present at a concentration of 1.5%, including:

- 0.5% in the rapid action digluconate form
- 1% in the prolonged action dihydrochloride form

PARAMETER OF EXCELLENCE

Chlorhexidine stands out as the primary reference element, both among natural substances and those synthesized in the laboratory, for its prominent role in controlling bacterial plaque and preventing gingivitis.

F. Graziani, S. Chieffi, C. Tomasi - Testo Atlante parodontologia e terapia implantare SIDP. implantare SIDP.

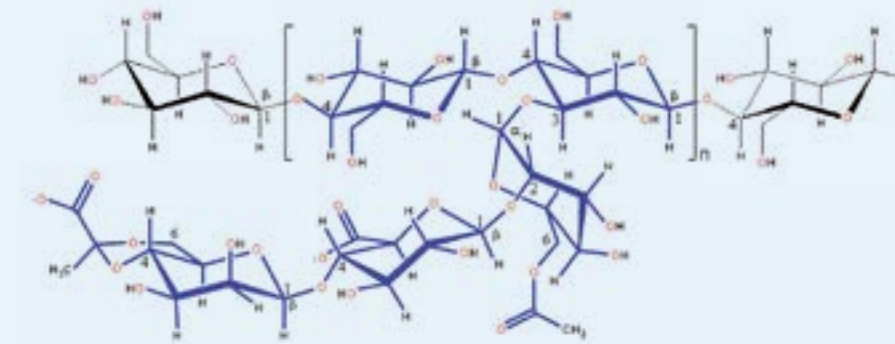

The recommended choice by the EUROPEAN FEDERATION OF PERIODONTOLOGY

Guidelines for the treatment of stage I - III periodontitis recommendations 2, 13, 14

The use of Chlorhexidine is the recommended choice for the treatment of stage I and III periodontitis, in addition to mechanical debridement. Chlorhexidine with prolonged release, when applied locally in patients with periodontitis, could be a solution to consider. No other chemical agent, except for some antibiotics, has been considered.

Sanz. M et al. Clinical Practice Guide Line - 2020 by EFP (Federazione Europea di Parodontologia).


Application of chlorhexidine pre and post oral surgery, periodontal surgery, and root planing procedures

The use of Chlorhexidine in the phases preceding and following surgical procedures offers the benefit of reducing bacterial presence in the oral cavity and preventing plaque formation in situations where mechanical hygiene may be difficult and inconvenient.

Newman e Addy.

Mucosal irritation - concentration-dependent idiosyncratic reaction - can only occur with the use of Chlorhexidine at high concentrations and for a prolonged period. Therefore, it can be concluded that Chlorhexidine at normal concentrations and usage methods cannot exert cytotoxic action.

WHY CHLO-SITE?

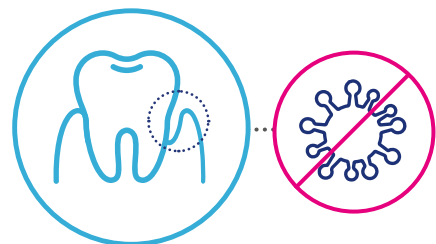

Unlike systemic pharmacological regimens, Chlo-Site can reach concentrations 100 times higher at the base of the pocket, **without the risk of toxicity and reducing the risk of developing tolerance or bacterial resistance.**

It is **READY TO APPLY**, suitable for use by dental hygienists as well.




The application time is less than **5 SECONDS.**



It remains within the pocket for **20+ DAYS**

It has a **FOLDABLE ATRAUMATIC** needle with dual lateral outlets.

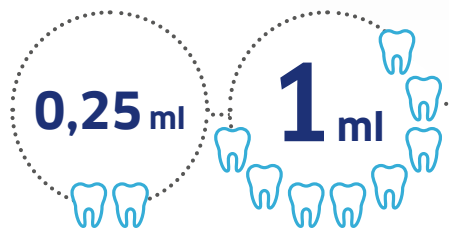

It is available in the **0,25 ml** version (single-use for half pockets) and the **1 ml** version for multiple pockets.

Bi-Chlorexidine 1.5%:

- **0,5% digluconate for rapid action**
- **1% in the prolonged action dihydrochloride form.**

Its mucoadhesive gel is **ADAPTABLE** to all types of pockets, even those smaller than 5 mm.



IT IS A CE MEDICAL DEVICE.




SCIENTIFIC EVIDENCE

Clinical, microbiological, and biochemical results of subgingival application of a **xanthan gel with chlorhexidine** in the treatment of periodontitis: a multicentric randomized study.

Randomized controlled split-mouth multicentric study.

Preliminary visit with 98 patients with moderate to advanced periodontitis.

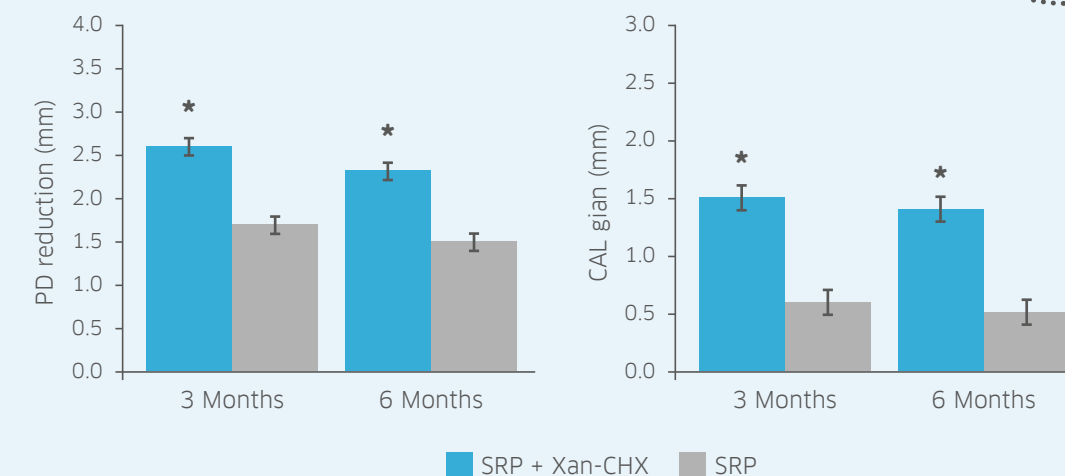

RESULTS:

The differences between the two treatment groups in PD reduction were 0.87 millimeters at 3 months and 0.83 millimeters at 6 months. Statistical analysis shows that the difference between the two groups is high and clinically relevant: P < 0.001.

The same considerations apply to CAL: 0.94 and 0.90 millimeters, respectively, at 3 and 6 months (P < 0.001).

It is worth noting that this improvement was also similar in the subgroup with pockets > 7 millimeters.

In other words, in cases with very deep pockets, **the combined use of CHLO-SITE with mechanical cleaning is capable of reducing periodontal disease to levels of mild to moderate severity.**

**Clinical and microbiologic effects of subgingival controlled-re-lease delivery of chlorhexidine chip in the treatment of periodontitis: a multicenter study. Paolantonio M, D'Angelo M, Grassi RF, Perinetti G, Picco-Iomini R, Pizzo G, Annunziata M, D'Archivio D, D'Ercole S, Nardi G, Guida L J Periodontol 2008; 79 (2): 271-82.



Therapeutic Periodontal Supplement
Your best ally for optimizing performance and costs

		CHLO-SITE	ANTIBIOTICS	CHIP	OTHER GEL
	GINGIVITIS Reduced Periodontium	✓ Supragingival Cleaning + CHLO-SITE	X	X	X
	GINGIVITIS Stable Periodontal Patient	✓ Supragingival Cleaning + CHLO-SITE	X	X	X
	PERIODONTITIS Stage I	✓ Supragingival Cleaning + CHLO-SITE	X	X	X
	PERIODONTITIS Stage II	✓ Supragingival Cleaning + CHLO-SITE	X	X	X
	PERIODONTITIS Stage III	✓ Subgingival Cleaning + CHLO-SITE	✓	✓	✓
	PERIODONTITIS Stage IV	✓ Subgingival Cleaning + CHLO-SITE	✓	✓	✓
	MUCOSITIS	✓ Element Cleaning + CHLO-SITE	X	X	✓
	PERI-IMPLANTITIS	✓ Subgingival Cleaning + CHLO-SITE	✓	✓	✓
	PROSTHETIC IMPLANTOLOGY	✓ CHLO-SITE	✓	X	✓
	PRE-PERIODONTAL SURGERY	✓ CHLO-SITE	✓	✓	✓



Unlike all other products
Chlo-Site is the only product useful in all periodontal issues.

Furthermore, thanks to the presence of chlorhexidine, it is the only gel that **guarantees effective.**

TREATMENTS FOR

Therapeutic Periodontal Supplement

The aforementioned plan should be considered as advice to the clinician based on current guidelines and the experience that our clinicians and users have had using Chlo-Site. It is not a protocol, and each clinician is free to use the product as they see fit based on their case-by-case evaluation.



1-2 sessions+

MUCOSITIS AND GINGIVITIS:

After site cleaning and thorough drying, apply Chlo-Site. In case of persistence, repeat the application in the next session. The decision to recall the patient and the interval between sessions depend on the clinician's assessment based on the PSR - Periodontal Screening and Recording Index.



3 sessions+

STAGE I PERIODONTITIS OR POCKETS LESS THAN 3 MM:

After removal of supragingival plaque, apply the product and recall the patient to repeat the treatment at least 2 times until the clinician deems it appropriate.



4 sessions+

STAGE II PERIODONTITIS OR 4 MM POCKETS:

After removal of subgingival plaque (SRP), apply the product and recall the patient to repeat the treatment at least 3 times until the clinician deems it appropriate.



5 sessions+

STAGE III PERIODONTITIS OR POCKETS GREATER THAN 5 MM:

After removal of subgingival plaque (SRP), apply the product and recall the patient to repeat the treatment at least 4 times until the clinician deems it appropriate.



6 sessions+

STAGE IV PERIODONTITIS OR POCKETS OVER 6 MM:

After removal of subgingival plaque (SRP), apply the product and recall the patient to repeat the treatment at least 5 times until the clinician deems it appropriate.



TREATMENT OF STABLE PERIODONTAL PATIENTS:

In patients with reduced periodontium and borderline probing situations, after removal of supragingival plaque, apply the product and recall the patient to repeat the treatment in all scheduled review sessions during the year.



PROSTHETIC IMPLANTOLOGY:

Use Chlo-Site in all implant cavities that need to be closed with prosthetic screws, both in the implant engagement (fixture) and on supra-implant prosthetic screws such as MUA (Multi Unit Abutment).



PRE-PERIODONTAL SURGERY:

In all situations where tissues need to be prepared for major periodontal or implant surgery, after supragingival or juxta-gingival debridement, apply the product into the sulcus or pocket to achieve better quality soft tissues.

* After applying Chlo-Site, ask the patient not to rinse for the next 2 hours.

With Chlo-Site you can:

- ✓ **SCHEDULE REGULAR CHECK-UPS** allowing the patient to have greater protection and loyalty.
- ✓ **INCREASE HEALING %** for all treatments, providing prolonged protection after professional cleaning.
- ✓ **SAVE** with packages that allow the use of Chlo-Site in all pockets, for every need.



Which package offer do you prefer?

6 Chlo-Site

5+1 STANDARD PACKAGE

- 5 packs of Chlo-Site 6 syringes 0.25 ml
- 1 pack of Chlo-Site 1 syringe 1 ml
- Total - 8,5 ml
- Number of pockets 68

7+2 GOLD PACKAGE

- 7 packs of Chlo-Site 6 syringes 0.25 ml
- 2 packs of Chlo-Site 1 syringe 1 ml
- Total - 12,5 ml
- Number of pockets 100

10+5 PREMIUM PACKAGE

- 10 packs of Chlo-Site 6 syringes 0.25 ml
- 5 packs of Chlo-Site 1 syringe 1 ml
- Total - 20 ml
- Number of pockets 160

4 Chlo-Site

5+1 STANDARD PACKAGE

- 5 packs of Chlo-Site 4 syringes 1 ml
- 1 pack of Chlo-Site 6 syringes 0.25 ml
- Total - 21 ml
- Number of pockets 172

7+2 GOLD PACKAGE

- 7 packs of Chlo-Site 4 syringes 1 ml
- 2 packs of Chlo-Site 6 syringes 0.25 ml
- Total - 31 ml
- Number of pockets 248

10+5 PREMIUM PACKAGE

- 10 packs of Chlo-Site 4 syringes 1 ml
- 5 packs of Chlo-Site 6 syringes 0.25 ml
- Total - 47,5 ml
- Number of pockets 380



BIBLIOGRAPHY

1. Evaluation of clinical efficacy therapeutic on the use of chlorhexidine gluconate to influence the shape of the gel and periochip during gingival inflammation in chronic periodontal disease. Saliu, Sahmedin and Kurtishi, Liburn and Nuhii, Nexhibe 2022 Sept.
2. Compare the Efficacy of Chlosite Gel as an Adjunctive Therapy after Scaling and Root Planing. Verma N, Saimbi CS, Gupta S, Kumar A, Tripathi AK. 2022 Jun.
3. Comparative Evaluation of Cost-Effective-ness, Clinical and Microbiological Parameters of Systemic Anti-biotics Versus Local Drug Delivery in Aggressive Periodontitis. Nandan B, Barman Roy D, Pant VA, Gupta V, Bhaduria U, Kaur H, Gupta O. Cureus. 2022 Jan.
4. Efficacy of Xanthan-Based Chlorhexidine Gel on the Levels of Interleukin-1 β in Chronic Periodontitis: An Interventional Study. Sajna HR, Ramesh A, Kedlaya MN, Thomas B. J Int Soc Prev Community Dent. 2021 Jul.
5. Treatment of periodontal disease using xanthan-based chlorhexidine gel. Abhishek Gautam, Kumar Manish, Raju Kumar. 2021 Feb.
6. Review on local drug delivery systems for periodontitis. K Kirshnananda Kamath, AR Shabaraya, M Kumar - 2021 Jan.
7. Professionally Delivered Local Antimicrobials in the Treatment of Patients with Periodontitis-A Narrative Review. Sholapurkar A, Sharma D, Glass B, Miller C, Nimmo A, Jennings E. Dent J (Basel). 2020 Dec 22;9(1):2.
8. Evaluation of different local drug delivery systems in the management of chronic periodontitis: a comparative study. Anna Abraham¹, Rekha Raghavan², Ajesh Joseph³, MP Shyamala Devi⁴, Megha Varghese⁵, 2020.
9. Adjunctive subgingival application of Chlorhexidine gel in nonsurgical periodontal treatment for chronic periodontitis: a systematic review and meta-analysis. Zhao, H., Hu, J. & Zhao, L. BMC Oral Health 20, 34 (2020).
10. Adjunctive effect of locally delivered antimicrobials in periodontitis therapy: A systematic review and meta-analysis. David Herrera, Paula Matesanz, Conchita Martín, Valerie Oud, Magda Feres, Wim Teughels.
11. Treatment of stage I-III periodontitis-The EFP S3 level clinical practice guideline. Sanz M, Herrera D, Kebschull M, Chapple I, Jepsen S, Be-glundh T, Sculean A, Tonetti MS; EFP Workshop Participants and Methodological Consultants. J Clin Periodontol. 2020 Jul;47 Suppl 22(Suppl 22):4-60.
12. Chlorhexidine gel used as antiseptic in periodontal pockets. Mummolo S, Severino M, Campanella V, Barlattani A Jr, Quinzi V, Marchetti E. J Biol Regul Homeost Agents. 2019 May-Jun;33(3 Suppl. 1):83-88. DENTAL SUPPLEMENT.
13. Chlorhexidine Gel Use in the Oral District: A Systematic Review. Luca Fiorillo Gels 2019, 5(2), 31; <https://doi.org/10.3390/gels5020031>.
14. Commercial local pharmacotherapeutics and adjunctive agents for nonsurgical treatment of periodontitis: A contemporary re-view of clinical efficacies and challenge. OL Tan, SH Safii, M Razali - Antibiotics, 2019 - mdpi.com.
15. Local drug delivery in the treatment of periodontal diseases. Ashika Sulthana, RT Arun, S Krishnaraj, Rajasekar Sundaram, 2019.
16. Therapeutic efficacy of chlorhexidine-gel on the depth of the periodontal pocket in chronic periodontal disease. Sali, Sahmedin; Fidoski, Jasmin Dec2018, Vol. 11 Issue 1, p92-95. 4p.
17. Physical key properties of antibiotic-free, PLGA/HPMC-based in-situ forming implants for local periodontitis treatment. K Agossa, M Lizambard, T Rongthong, E Delcourt-Debruyne, J Siepmann, F Siepmann, 2017 Apr 15;521(1-2):282-293. doi: 10.1016/j.ijpharm.2017.02.039. Epub 2017 Feb 20.
18. The clinical and metabolic effects of subgingival application of xanthan-based chlorhexidine gel in Type 2 diabetic patients with chronic periodontitis. Faramarzi M, Shirmohammadi A, Chitsazi M, Shamami MS, Ghanitab S. Dent Res J (Isfahan). 2017 Sep-Oct;14(5):299-305. doi: 10.4103/1735-3327.215961. PMID: 29109743; PMCID: PMC5654223.
19. Evaluation of clinical efficacy of 0.2% chlorhexidine irrigation, 1.5% chlorhexidine gel and 2.5mg biodegradable chlorhexidine chip as an adjunct to scaling and root planing in the management of Chronic Periodontitis. Shubhra Vaish, Vidya Dodwad, Aakriti Mahajan, Sonakshi Gupta
20. Comparative evaluation of subgingivally delivered xanthan-based chlorhexidine gel and herbal extract gel in the treatment of chronic periodontitis. Phogat M, Rana T, Prasad N, Baiju CS. J Indian Soc Periodontol. 2014 Mar;18(2):172-7. doi: 10.4103/0972-124X.131319.
21. Oral antiseptic and periodontitis: a clinical and microbiological study. Mummolo S, D'Ercole S, Marchetti E, Campanella V, Martinelli D, Marzo G, Tripodi D. Oral Health Dent Manag. 2014 Sep;13(3):698-702. PMID: 25284540.
22. Efficacy of xanthan-based chlorhexidine gel as an adjunct to scaling and root planing in treatment of the chronic periodontitis. Meetu J, Deepak D, Preet J, Balaji M, Bipin Y, Neema S. J Indian Soc Periodontol. 2013; 17 (4): 439-43.
23. Efficacy of Subgingival Application of Xanthan-based Chlorhexidine Gel Adjunctive to Full-mouth Root Planing Assessed by Real-time PCR: A Microbiologic and Clinical Study. Chitsazi MT, Kashefimehr A, Pourabbas R, Shirmohammadi A, Ghasemi Barghi V, Daghigh Azar B. J Dent Res Dent Clin Dent Prospects. 2013;7(2):95-101. doi: 10.5681/joddd.2013.017. Epub 2013 May 30. PMID: 23875088; PMCID: PMC3713868.
24. Effect of insertion of xanthan-based chlorhexidine gel in the maintenance phase following the treatment of chronic periodontitis. Verma A, Sanghi S, Grover D, Aggarwal S, Gupta R, Pandit N J Indian Soc Periodontol. 2012; 16 (3): 381-5.
25. A randomized clinical trial on the clinical and microbiological efficacy of a xanthan gel with chlorhexidine for subgingival use. Matesanz P, Herrera D, Echeverría A, O'Connor A, González I, Sanz M. Clin Oral Invest - Springer-Verlag 2012.
26. The comparative evaluation of xanthan gel with chlorhexidine (CHLO-SITE) in smokers and non-smokers: a clinical and microbiological assessment. Chandra C, Narayan Valavalkar N, Vandana KL. J Indian Soc Periodontol 2011; 15 (3).
27. Effect of CHLO-SITE (xanthan gel with chlorhexidine) on clinical and microbiological parameters in smokers - a case series. Chandra C, Chandra S. Indian Journal of Forensic Medicine and Toxicology. 2010; 4 (10): 370-3.
28. Management parodontale mediante antisettici a rilascio controllato di clorexidina. Mummolo S, Marchetti E, Di Martino S, D'Aversa L, Marzo G. DoctorOs 2010; 21 (9): 1139-43.
29. Clinical, microbiological and biochemical effects of subgingival administration of a xanthan-based chlorhexidine gel in the treatment of periodontitis: a randomized multicenter trial. Paolantonio M, D'Ercole S, Pilloni A, D'Archivio D, Lisanti L, Graziani F, Femminella B, Sammartino G, Perillo L, Tetè S, Perfetti G, Spoto G, Piccolomini R, Perinetti G. J Periodontol 2009; 80: 1479-92.
30. Comparative evaluation of subgingival delivered 10% doxycycline hyclate and xanthan-based chlorhexidine gels in the treatment of chronic periodontitis. Gupta R, Pandit N, Aggarwal S, Verma A. J Contemp Dent Pract 2008; (9) 7: 25-32.
31. Clinical and microbiologic effects of subgingival controlled-release delivery of chlorhexidine chip in the treatment of periodontitis: a multicenter study. Paolantonio M, D'Angelo M, Grassi RF, Perinetti G, Picco-lomini R, Pizzo G, Annunziata M, D'Archivio D, D'Ercole S, Nardi G, Guida L J Periodontol 2008; 79 (2): 271-82.
32. Internal decontamination of dental implants: an in vivo randomized microbiologic 6-month trial on the effects of a chlorhexidine gel. Paolantonio M, Perinetti G, D'Ercole S, Graziani F, Ca-tamo G, Sammartino G, Pic-colomini R. J Periodontol 2008; 79: 1419-25.
33. Gingival crevicular fluid alkaline phosphatase activity reflects periodontal healing/re-current inflammation phases in chronic periodontitis patients. Perinetti G, Paolantonio M, Femminella B, Serra E, Spoto G. J Periodontol 2008; 79: 1200-7.
34. The effect of locally delivered xanthan-based CHLO-SITE gel with scaling and root planning in the treatment of chronic periodontitis: microbial findings. Abrishami M, B Iramloo B, Ansari G, Eslami G, Bagheban AA, Anaraki M. Dent Res J 2008; 5 (2): 47-52.
35. Non-surgical periodontal therapy using a novel chlorhexidine based xanthan-gel: a split-mouth study. Rusu D, Stratul SI, Necker A, Benta A. Int Poster J Dent Oral Med 2005; 7 (3): Poster 286.



Available in 3 formats

- ✓ **BOX OF 6 SYRINGES OF 0.25 ML**
- ✓ **BOX OF 4 SYRINGES OF 1 ML**
- ✓ **BOX OF 1 SYRINGE OF 1 ML**

Reserved for the professional body.



Want to learn more?
Contact us for further information.
+39 051 57 53 53
info@ghimas.it

GHIMAS S.p.A.

Via Cimarosa, 85 - 40033 Casalecchio di Reno (BO) - Italy
+39 051 57 53 53 - www.ghimas.com - info@ghimas.it

